



**Model: 507**

**Sonic Electric Toothbrush**

---

## **Instruction Manual**

Thanks for choosing 507. Before use, please carefully read the instructions and keep it at a safe place for future reference.

### **1. Features**

- Sonic vibration, over 40,000 strokes per minute, dynamic cleaning.
- Remove up to 100% more plaque within one week. \*
- Improve gum health in just 2 weeks. \*
- Achieve noticeably whiter teeth in just 2 weeks. \*

- 5 modes: White, Clean, Sensitive, Polish and Massage.
- Unique design of tooth-waved bristles for thorough cleaning of gum and hard to reach areas.
- 2 minutes smart timer with interval pause reminding each 30 seconds.
- Replaceable toothbrush head.
- Waterproof rated IPX7.
- USB charging, easy and convenient for use at home or traveling.
- Long lasting battery, 1 charging for around 30 days use.  
\*compared with a regular manual toothbrush.

## 2. Designs

- ① Toothbrush head
- ② "Power" button
- ③ White
- ④ Clean
- ⑤ Sensitive
- ⑥ Polish
- ⑦ Massage
- ⑧ Charging light
- ⑨ Charging port



## 3. Operations


### ☆ Using the toothbrush

- Attach the brush head and make sure it is locked in place.
- Choose your preferred mode by pressing "Power" button. You may shift among the 5 modes "White, Clean, Sensitive, Polish and Massage" by continuous pressing.
- Wet the bristles and apply appropriate amount of toothpaste.
- Place the bristles with toothpaste against the teeth at a slight angle towards the gumline.
- Press the "Power" button to turn on.
- Gently move the brush head across the teeth on front, back and chewing surfaces. Do not apply too much pressure to the teeth or apply the bristles to one area for long time.
- 2 minutes is recommended by dentists to brush the teeth and there is interval pause reminding after each 30 seconds.
- Move to next section when you feel the pause and the toothbrush will automatically power off at 2 minutes.
- To avoid splashing toothpaste, only turn on the toothbrush after you have placed it against your teeth and only remove it from the mouth after it is turned off.


### ☆ Changing the toothbrush head

- Install: Slide toothbrush head down to the handle. Turn clockwise until A meets T.
- Detach: Turn toothbrush head counter clockwise until A meets B. Pull it out along the handle.

Install

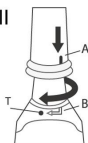

Detach

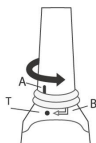

## ☆Charging the toothbrush

- When battery power level is low, the charging indicator light will flash 10 times, to remind charging.
- Connect the toothbrush with USB cable to any power source to start charging with the charging indicator light on.
- When the battery is fully charged, indicator light will turn off automatically and the charging unit will shut off automatically to avoid overcharging.
- For the first time use or when it is not used for long time, it is normal that the indicator does not light for initial 5 minutes.

## 4. Specifications

Toothbrush power source: DC3.7V

Product material: ABS/TPE

Bristle: Nylon

Toothbrush power consumption: 1.0W

Working temperature: 0-50℃

## 5. Precautions

- Do not disassemble, repair, or alter the product without authorization.
- People who have severe periodontal disease or oral surgery should consult a doctor before use.
- Toddlers, infants and people who are unable to operate should be supervised by a parent or adult.
- Do not soak the product in water for long time.
- Do not place the toothbrush in high temperature condition or expose it to the sun.
- Do not force the bristles hard against the teeth or gums


## 6. After-sales service

- 1-year warranty is offered, since the day of purchase. You may contact Fairywill.
  - During warranty period, for any defects arisen from materials or product process, we will replace the parts or the whole product according to specific circumstances.
  - When you apply for free repair, you must submit 507 and proof which can prove the date of purchase (purchasing invoice) to Fairywill.
  - The following conditions are not included in the warranty:
    - Damages caused by improper use, for example, wrong use of power supply, plugged into improper socket, artificial damage etc.
    - Ageing and wearing-out caused by normal use.
    - Problems occur due to its use for other purposes, or improper use or use of components not supplied by Fairywill.
-

# Fairywill

## Smart Oral Irrigator (F5020E) Instruction Manual

Thanks for choosing Fairywill Oral irrigator. Before operating the product, please read the manual carefully and keep it for future reference.

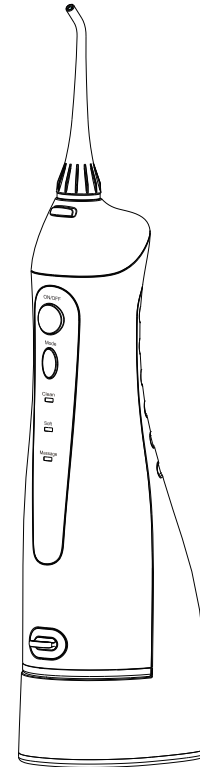

ISO9001 & ISO13485 Certificated Manufacturer



Printed in China



## 1. Important Safeguards:

All information in this section is provided to instruct you on the correct and safe operation of this device so that you can avoid injury to yourself and others as well as damage to property.


- Toddlers, infants, people who are unable to operate the device should not use it by themselves. Please also keep it away and out of reach from children.
- Never try to alter, disassemble or repair the device.
- Do not place or store the product where it will fall or end up in a tub or sink.
- Do not reach for the product if it has fallen into water. Unplug immediately.
- Do not use the product while bathing.

## 2. Notes:

- This appliance can be used by children from age 8 and above; and people with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the device in a safe way and understand the hazards involved.
- Children shall not play with the device. Cleaning and user maintenance shall not be made by children without supervision.
- The device can be only be used with the power supply unit provided with the device.
- The battery in the device must be removed from the device before it is discarded.
- The device must be disconnected from the powersupply when removing the battery.
- The battery should be disposed of safely.
- Regarding the removal of the battery, please see "Disposal" on page 7.
- Do not fill the tank with water warmer than 40 °C.
- If the external flexible cable or cord of this transformer is damaged, it shall be replaced by the manufacturer, a service agent or a similarly qualified person in order to avoid any Hazard.
- Keep the cord away from heated surfaces.
- Please do not use it if you lose the nozzle.
- Do not submerge the device in water.
- Do not use for any purpose other than oral cleaning.
- Keep nozzle away from young children. The nozzle is small enough to swallow and can be lodged in the throat.

- It is not recommended to use mouthwash as it is corrosive and may affect the lifespan of the irrigator.
- People who experience tooth or gum pain may not be able to use this device and should consult with their dentist before use.
- Do not press the nozzle too firmly against your teeth or gums, it will hurt your gums if you do so.

## 3. Designs:

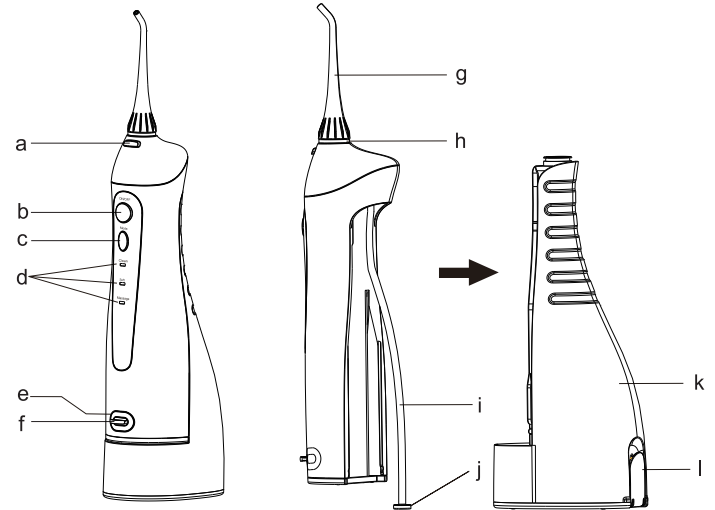

- a. Nozzle release button
- b. On/Off button
- c. Mode button
- d. Mode indicator
- e. Charging base
- f. Charging base cover

- g. Nozzle
- h. Identification ring
- i. Nozzle suction tube
- j. Filter
- k. Water tank (300ML)
- l. Water tank cap

### Note:

There may be some liquid remaining inside the main device before use (tank, water suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

## 4. How To Use:

### 4.1 Charging the device ( As Fig. 1)

- Please cover the lid of the tank tightly before using the water flosser to avoid water leakage.
- Take the charger base cover down and insert the device socket (make sure no liquid is remaining inside the device plug).
- Insert the charging adaptor securely into the power outlet.
- The “CH” indicator blinks when the device is charging, and it will turn off after the device is full charged.
- Put the cover on the charger base again after charging is finished.

#### Note:

If the “CH” indicator is blinking during or after an irrigation, that means the battery needs to be recharged. A full charge will take about 4 hours.

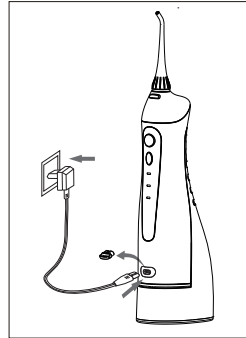


Fig. 1

### 4.2 Use the Oral Irrigator

- Attach the nozzle onto the unit. (Fig. 2)
  - Open the cap of the water tank to fill the tank with water which should not exceed 40°C, then close the cap.(Fig. 3)
  - Or, slide the water tank down off the unit, fill the tank with water which should not exceed 40°C, then attach the tank on the unit. (Fig. 4)
  - Press the “Mode” key to select the mode you desired. During irrigating, you can change the working mode. (Fig. 5)
- Clean:** for normal cleaning, it will last for around 90s (+/-15s).
- Soft:** for gentle cleaning of sensitive gums for around 110s(+/-15s).
- Massage:** for cleaning and massaging the gums for around 100s(+/-15s).
- Hold the unit in an upright position with the nozzle tip pointing toward your teeth.
  - Keep your mouth opened slightly so that the water can flow out.

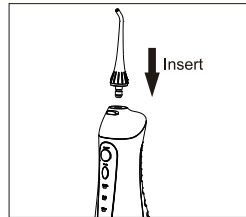


Fig. 2

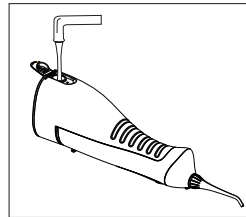


Fig. 3

- Press “On/Off” key to start irrigating. (Fig. 6)
  - 1) Direct the water jet flow at a 90-degree angle to the teeth and gums.
  - 2) Direct the water jet slowly along the teeth.
  - 3) Direct the water jet along the gum line with the nozzle tip at a 90 degree angle to the gum tissue.
  - 4) Cleaning between the teeth and around any braces, crowns or bridgework thoroughly.
- Press the “On/Off” key and then press the nozzle release button to replace other nozzles when needed. (Fig.7)
- Once irrigating is complete, press the “On/Off” key again to turn off the device.

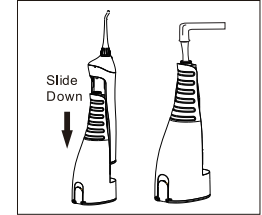


Fig. 4

### 4.3 After using the appliance (Fig.8)

- Open the cap of the water tank and dispose of the water in the tank.
- Press the “On/Off” key to turn on the unit, and drain the water inside the device.
- Press the “On/Off” key again to turn off the unit, and close the water tank cap.
- Wipe the device with a dry cloth.

#### Note:

To prevent bacteria from propagating, make sure there is no remaining water inside the unit. If you're not going to use this appliance for a long time, be sure to wipe and dry the appliance as soon as possible before storing it.

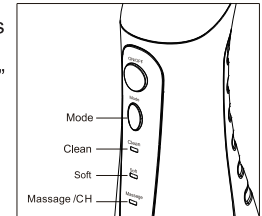


Fig. 5

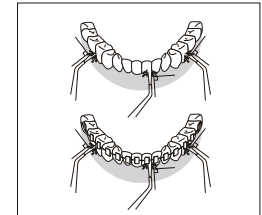


Fig. 6

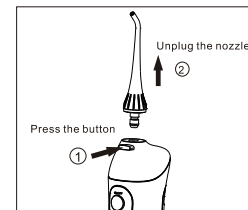


Fig. 7

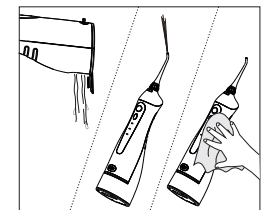


Fig. 8

## 5. Maintenance:

The device can be only cleaned by water or a neutral cleaning agent. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not use hot water over 50°C to wash it. Keep the device away from high-temperature areas and avoid direct sunlight.

### 5.1 Main Unit

- Remove the dirt by using a cloth with liquid soap to clean the unit and a dry cloth to wipe.
- When cleaning the main unit, do not submerge it in water.
- Dry the unit after washing with water.

### 5.2 Water Tank Cap

Wipe off any grime and dirt with a cloth.

### 5.3 Nozzle

- Wash with water and wipe with a soft cloth.
- Replace the nozzle every 6 months.

### 5.4 Water Suction Tube

- Wash with water and wipe with a soft cloth.
- Do not bend, pull or twist the tube.

### 5.5 Tank

- Wash with water and wipe with soft cloth.
- Dry the tank when it is not in use for an extended period of time.


## 6. Troubleshooting:

Problem	Possible Cause	Solution
It does not work	You've just purchased the product or it has been left unused for more than 3 months.	Charge the unit for at least 4 hours.
It only works for a few minutes even after charging	The battery has reached the end of its lift cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the unit for at least 4 hours.
Water pressure is very low	There is dirt stuck in the nozzle	Replace the nozzle
	The nozzle is deformed	
	There is dirt stuck in the filter	Clean the filter
Water does not come out	The tank is empty	Fill up tank
	The appliance is tilted too much while in use	Use with the appliance in a vertical state

## 7. Disposal:

The rechargeable battery inside your appliance is easy to remove for recycling. At the end of the product's life and before disposal, please remove the rechargeable battery from the handle by following the instructions below. During the operation below, be sure to protect your eyes, hands, fingers, and the surface on which you work.

- Remove 2 pieces of rubber feet on the bottom of the handle (Like ① on Fig. 9)
- Insert a screwdriver into the slot located at the bottom of the handle and turn counter-clockwise to release the bottom cap. (Like ② on Fig. 9)
- Slide the water tank down off the device and unscrew the two screws next to the suction tube. (Fig. 10)

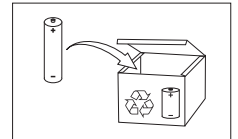

- Open the top cover and put the inner part outside of the unit, then unscrew all the screws of the inner part.(Fig. 11)
- Cut off the battery line and pull out the batteries from the inner part.(Fig. 12)

If you have recycling questions, please contact your local waste management office.

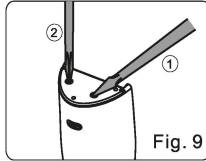


Fig. 9

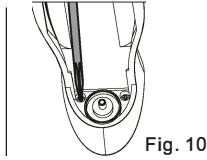


Fig. 10

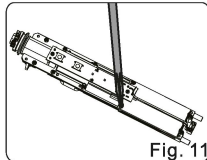


Fig. 11

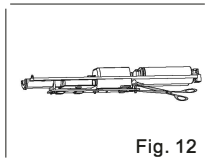


Fig. 12

## 8. Warranty:

- 1-year warranty is offered, starting from the first day of purchase. You may contact Fairywill.
- During the warranty period, for any defect that may arise from materials or product process, we will replace the parts or the whole product according to certain circumstances.
- When you apply for free repair, you must submit the product and proof that can prove the date of purchase to Fairywill.
- The following conditions are not included in the warranty:
  - Damages caused by improper use, for example, wrong use of power supply, plugged into an improper socket, artificial damage etc.
  - Ageing and normal wear and tear from use.
  - Problems occur due to its use for other purposes, or improper use or use components not supplied by Fairywill.
- This warranty does not affect your statutory rights.

## 9. Contact:

- 📍 Fairywill
- 🌐 [www.fairywill.com](http://www.fairywill.com)
- ✉ [support@fairywill.com](mailto:support@fairywill.com)

Made in China

Distributed by:

Zhuhai Zonghengyoupin Kuajingdianshang  
Youxiangongsi, Zhuhai,  
Guangdong, 519031, CN



ISO9001 & ISO13485 Certified Manufacturer



Printed in China