



**Model: 507**

**Sonic Electric Toothbrush**

---

## **Instruction Manual**

Thanks for choosing 507. Before use, please carefully read the instructions and keep it at a safe place for future reference.

### **1. Features**

- Sonic vibration, over 40,000 strokes per minute, dynamic cleaning.
- Remove up to 100% more plaque within one week. \*
- Improve gum health in just 2 weeks. \*
- Achieve noticeably whiter teeth in just 2 weeks. \*

- 5 modes: White, Clean, Sensitive, Polish and Massage.
- Unique design of tooth-waved bristles for thorough cleaning of gum and hard to reach areas.
- 2 minutes smart timer with interval pause reminding each 30 seconds.
- Replaceable toothbrush head.
- Waterproof rated IPX7.
- USB charging, easy and convenient for use at home or traveling.
- Long lasting battery, 1 charging for around 30 days use.  
\*compared with a regular manual toothbrush.

## 2. Designs

- ① Toothbrush head
- ② "Power" button
- ③ White
- ④ Clean
- ⑤ Sensitive
- ⑥ Polish
- ⑦ Massage
- ⑧ Charging light
- ⑨ Charging port



## 3. Operations


### ☆ Using the toothbrush

- Attach the brush head and make sure it is locked in place.
- Choose your preferred mode by pressing "Power" button. You may shift among the 5 modes "White, Clean, Sensitive, Polish and Massage" by continuous pressing.
- Wet the bristles and apply appropriate amount of toothpaste.
- Place the bristles with toothpaste against the teeth at a slight angle towards the gumline.
- Press the "Power" button to turn on.
- Gently move the brush head across the teeth on front, back and chewing surfaces. Do not apply too much pressure to the teeth or apply the bristles to one area for long time.
- 2 minutes is recommended by dentists to brush the teeth and there is interval pause reminding after each 30 seconds.
- Move to next section when you feel the pause and the toothbrush will automatically power off at 2 minutes.
- To avoid splashing toothpaste, only turn on the toothbrush after you have placed it against your teeth and only remove it from the mouth after it is turned off.


### ☆ Changing the toothbrush head

- Install: Slide toothbrush head down to the handle. Turn clockwise until A meets T.
- Detach: Turn toothbrush head counter clockwise until A meets B. Pull it out along the handle.

Install

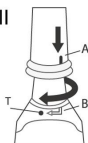

Detach

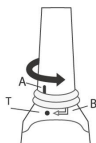

## ☆Charging the toothbrush

- When battery power level is low, the charging indicator light will flash 10 times, to remind charging.
- Connect the toothbrush with USB cable to any power source to start charging with the charging indicator light on.
- When the battery is fully charged, indicator light will turn off automatically and the charging unit will shut off automatically to avoid overcharging.
- For the first time use or when it is not used for long time, it is normal that the indicator does not light for initial 5 minutes.

## 4. Specifications

Toothbrush power source: DC3.7V

Product material: ABS/TPE

Bristle: Nylon

Toothbrush power consumption: 1.0W

Working temperature: 0-50℃

## 5. Precautions

- Do not disassemble, repair, or alter the product without authorization.
- People who have severe periodontal disease or oral surgery should consult a doctor before use.
- Toddlers, infants and people who are unable to operate should be supervised by a parent or adult.
- Do not soak the product in water for long time.
- Do not place the toothbrush in high temperature condition or expose it to the sun.
- Do not force the bristles hard against the teeth or gums

## 6. After-sales service

- 1-year warranty is offered, since the day of purchase. You may contact Fairywill.
  - During warranty period, for any defects arisen from materials or product process, we will replace the parts or the whole product according to specific circumstances.
  - When you apply for free repair, you must submit 507 and proof which can prove the date of purchase (purchasing invoice) to Fairywill.
  - The following conditions are not included in the warranty:
    - Damages caused by improper use, for example, wrong use of power supply, plugged into improper socket, artificial damage etc.
    - Ageing and wearing-out caused by normal use.
    - Problems occur due to its use for other purposes, or improper use or use of components not supplied by Fairywill.
-