



**Model: 507**

**Sonic Electric Toothbrush**

---

## **Instruction Manual**

Thanks for choosing 507. Before use, please carefully read the instructions and keep it at a safe place for future reference.

### **1. Features**

- 
- Sonic vibration, over 40,000 strokes per minute, dynamic cleaning.
  - Remove up to 100% more plaque within one week. \*
  - Improve gum health in just 2 weeks. \*
  - Achieve noticeably whiter teeth in just 2 weeks. \*

- 5 modes: White, Clean, Sensitive, Polish and Massage.
- Unique design of tooth-waved bristles for thorough cleaning of gum and hard to reach areas.
- 2 minutes smart timer with interval pause reminding each 30 seconds.
- Replaceable toothbrush head.
- Waterproof rated IPX7.
- USB charging, easy and convenient for use at home or traveling.
- Long lasting battery, 1 charging for around 30 days use.  
\*compared with a regular manual toothbrush.

## 2. Designs

- ① Toothbrush head
- ② "Power" button
- ③ White
- ④ Clean
- ⑤ Sensitive
- ⑥ Polish
- ⑦ Massage
- ⑧ Charging light
- ⑨ Charging port



## 3. Operations


### ☆ Using the toothbrush

- Attach the brush head and make sure it is locked in place.
- Choose your preferred mode by pressing "Power" button. You may shift among the 5 modes "White, Clean, Sensitive, Polish and Massage" by continuous pressing.
- Wet the bristles and apply appropriate amount of toothpaste.
- Place the bristles with toothpaste against the teeth at a slight angle towards the gumline.
- Press the "Power" button to turn on.
- Gently move the brush head across the teeth on front, back and chewing surfaces. Do not apply too much pressure to the teeth or apply the bristles to one area for long time.
- 2 minutes is recommended by dentists to brush the teeth and there is interval pause reminding after each 30 seconds.
- Move to next section when you feel the pause and the toothbrush will automatically power off at 2 minutes.
- To avoid splashing toothpaste, only turn on the toothbrush after you have placed it against your teeth and only remove it from the mouth after it is turned off.


### ☆ Changing the toothbrush head

- Install: Slide toothbrush head down to the handle. Turn clockwise until A meets T.
- Detach: Turn toothbrush head counter clockwise until A meets B. Pull it out along the handle.

Install

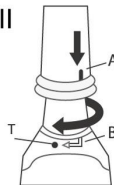

Detach

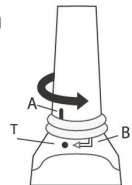

### ☆Charging the toothbrush

- When battery power level is low, the charging indicator light will flash 10 times, to remind charging.
- Connect the toothbrush with USB cable to any power source to start charging with the charging indicator light on.
- When the battery is fully charged, indicator light will turn off automatically and the charging unit will shut off automatically to avoid overcharging.
- For the first time use or when it is not used for long time, it is normal that the indicator does not light for initial 5 minutes.

### 4. Specifications

Toothbrush power source: DC3.7V

Product material: ABS/TPE

Bristle: Nylon

Toothbrush power consumption: 1.0W

Working temperature: 0-50°C

### 5. Precautions

- Do not disassemble, repair, or alter the product without authorization.
- People who have severe periodontal disease or oral surgery should consult a doctor before use.
- Toddlers, infants and people who are unable to operate should be supervised by a parent or adult.
- Do not soak the product in water for long time.
- Do not place the toothbrush in high temperature condition or expose it to the sun.
- Do not force the bristles hard against the teeth or gums

# Fairywill

**Model: 2001**  
**Sonic Electric Toothbrush**




Thank you for choosing Fairywill 2001 sonic toothbrush. Please read the instructions carefully before using the product, and keep the user manual in a safe place for future reference.

## 1. Features:


- Cleans your teeth and gums effectively with a vibration frequency of 38,000 times per minute.
- Includes 3 modes (low, high and massage) to meet different demands.
- Anti-slipping – easy, firm grip and ergonomically designed.
- Waterproof rating at IPX7.
- High-level safety as it works with durable Ni-MH battery. Battery life can last for 21 days on a full charge and has a low-battery indicator. A USB charging cable is also included with the product.
- A smart timer will remind you every 30 seconds to change your brushing area.
- The brush will shut off automatically after 2 minutes to keep an accurate brushing time.

## 2. Designs:


### Component Name

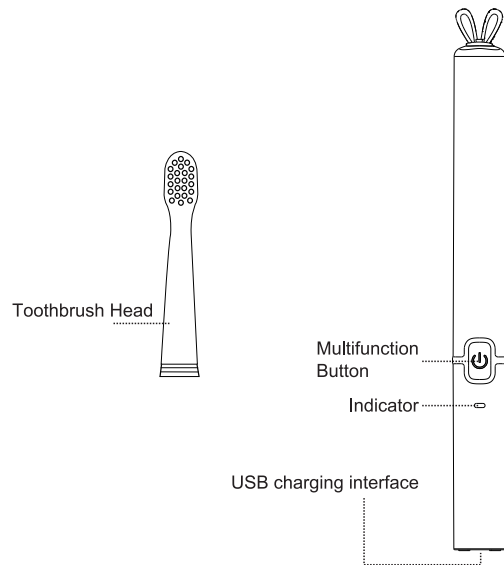

### Instruction

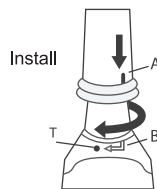

Install  
Slide toothbrush head down as direction A, turn clockwise let A meet T

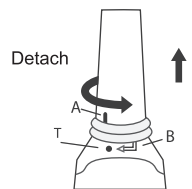

Detach  
Turn toothbrush head counter clockwise until A meet B, then pull out as direction

### **3. Important Safeguards:**

- When using the toothbrush for the first time, some people may need time to adapt to the brush and might experience feelings such as numbness. This is normal for first-time users. The feelings should fade away after several uses, and you should get used to it within a week.
- Do not open, repair or alter the product without seeking professional help.
- Do not hit or drop the product. If it is defective, please stop using it.
- Please check if the product is compatible with the local voltage in your area before charging the product for use.
- Only certain toothbrush heads will fit this product.
- Please consult your dentist before using it if you've had surgery in your gums or oral in the last two months.
- If you experience bleeding after using the toothbrush, or if the bleeding continues for more than one week, please consult your dentist immediately.
- The product is designed only for cleaning your teeth, gums and tongue, and not for other purpose.
- Children should be supervised when using the toothbrush to ensure proper and safe use. Please keep it away from children when not in use.
- Do not expose the toothbrush to high heat or directly under sunlight.

### **4. Specifications:**

Working temperature: 0-50 °C

Working voltage: DC1.2V

Bristles material: Nylon

Input voltage: DC5.0V

Input power: 0.6W

Working power: 0.55W

Charging time: ≤6 hours

### **5. How to use:**

- Connect the product to power supply by using the USB cable and charge the product.
- Attach the toothbrush head to the handle, and make sure it is locked in the correct place. Damp the bristles and apply toothpaste.
- Gently place the brush head on your teeth.
- Gently move the brush head across the teeth, brushing the front, back and chewing surfaces for 10 seconds each.
- Press the Multifunctional button to turn on the toothbrush and brush all surfaces of the teeth
- Press the Multifunctional button to change brushing modes:
  - Low: For gentler cleaning of your teeth, tongue and sensitive areas. (Indicator light in green)
  - High: For efficient cleaning of your teeth. (Indicator light in red)
  - Massage: For gentle stimulation, designed to keep your gums firm and healthy. (Indicator light in blue)

- Move the brush to the next area when the auto reminder pauses briefly every 30 seconds. After finishing all 4 areas, the brush will automatically shut off in 2 minutes.

## **6. Charging Instructions:**


- When the battery is low, the charging indicator light will flash 10 times.
- The indicator light will remain red during charging. Once fully charged, the light will turn green indicating the battery is full. Once the toothbrush is fully charged, the charging unit will shut off to prevent overcharging.
- For first-time usage or if the battery has been on low battery for a long period of time, it is normal if the indicator light does not flash in the first 5 minutes.

## **7. After-sales service:**

- 1-year warranty starts on the day of purchase, please contact Fairywill when needed.
- Within the warranty period, any defects arisen from materials or product process, we will replace the parts or the whole product depending on the condition.
- When you apply for a free repair, you must submit 2001 and proof of the date of purchase (purchasing invoice) to Fairywill.
- The following conditions are not included in the warranty:
  - Damages caused by improper use, for example, wrong use of power supply, plugged into improper socket, artificial damage, etc.
  - Ageing and wearing-out caused by normal use.
  - Damages due to product being used for other purposes, or use of components not supplied by Fairywill.

## **8. Contact:**

 Fairywill

 [www.fairywill.com](http://www.fairywill.com)

 [support@fairywill.com](mailto:support@fairywill.com)

Made in China

Distributed by:

Zhuhai Zonghengyoupin Kuajingdianshang Youxiangongsi, Zhuhai, Guangdong, 519031, CN